


Aroclor 1260: DIN 51527 on ZB-XLB-HT**Column:** Zebron[™] ZB-XLB-HT, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G024-11**Oven Profile:** 50 °C for 0.5 min to 210 °C @ 40 °C/min for 3 min to 230 °C @ 30 °C/min for 5 min to 250 °C @ 30 °C/min for 5 min to 320 °C @ 40**Carrier Gas:** Constant Flow Helium, 1.5 mL/min**Injection:** Split 2:1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (350°C)**Detector Info:** Pressure pulse used at 40 psi for first 0.25 min**Products used in this application:****Analyst Note:** Total concentration of Aroclors was 90 ppm in iso-octane

Aroclor 1260: DIN 51527 on ZB-XLB-HT

ANALYTES:

- 1 TCMX
- 2 BZ# 31
- 3 BZ# 28
- 4 BZ# 52
- 5 BZ# 101
- 6 BZ# 118
- 7 BZ# 153
- 8 BZ# 105
- 9 BZ# 138
- 10 BZ# 167
- 11 BZ# 156
- 12 BZ# 180
- 13 BZ# 170
- 14 BZ# 189
- 15 DCB

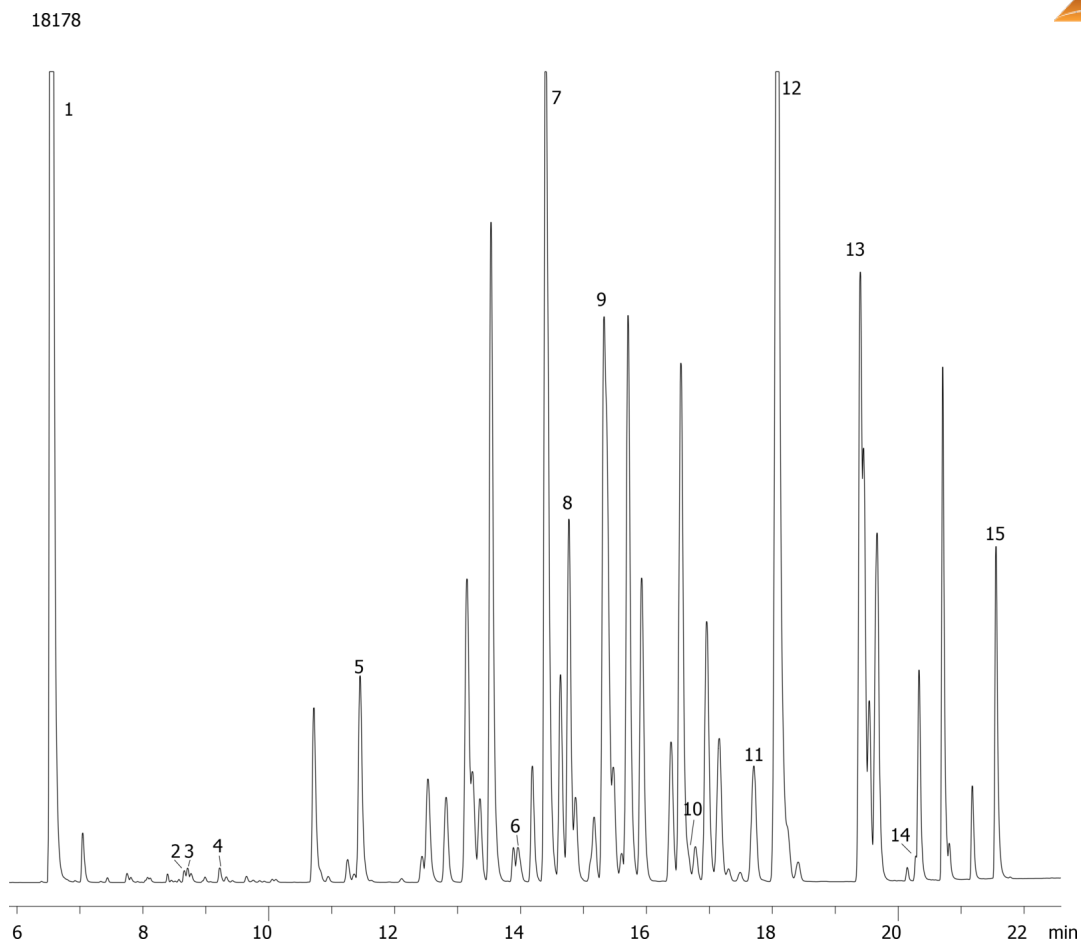