

Aflatoxins on Prodigy 5 ODS (2)

Column: Prodigy™ 5 µm ODS-2 150 Å, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-3300-E0
Elution Type: Isocratic
Eluent A: Water
Eluent B: CH₃OH
Gradient Profile:


Step No.	Time (min)	Pct A	Pct B
1	0	55	45

Flow Rate: 0.8 mL/min
Col. Temp.: ambient
Detection: Refractive Index @ 0.0000000000 (ambient)



Products used in this application:

Prodigy™

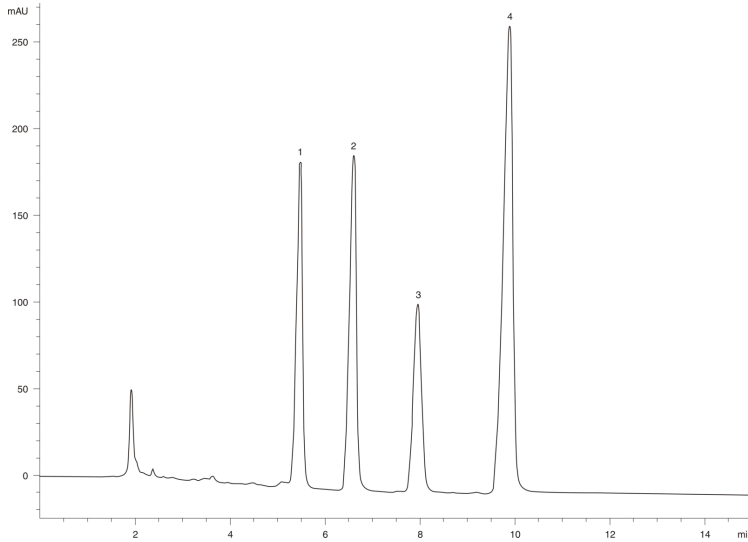

ANALYTES:

- 1 Aflatoxin G2
- 2 Aflatoxin G1
- 3 Aflatoxin B2
- 4 Aflatoxin B1

