

Felbamate

Column: Prodigy™ 5 µm ODS-3 100 Å, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4097-E0
Elution Type: Isocratic
Eluent A: 8,4 g KH₂PO₄ in 1700 ml water, pH 2.5 w/ H₃PO₄ + 300 ml acetonitrile
Gradient

Step No.	Time (min)	Pct A
1	0	100


Flow Rate: 2 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287
- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy™

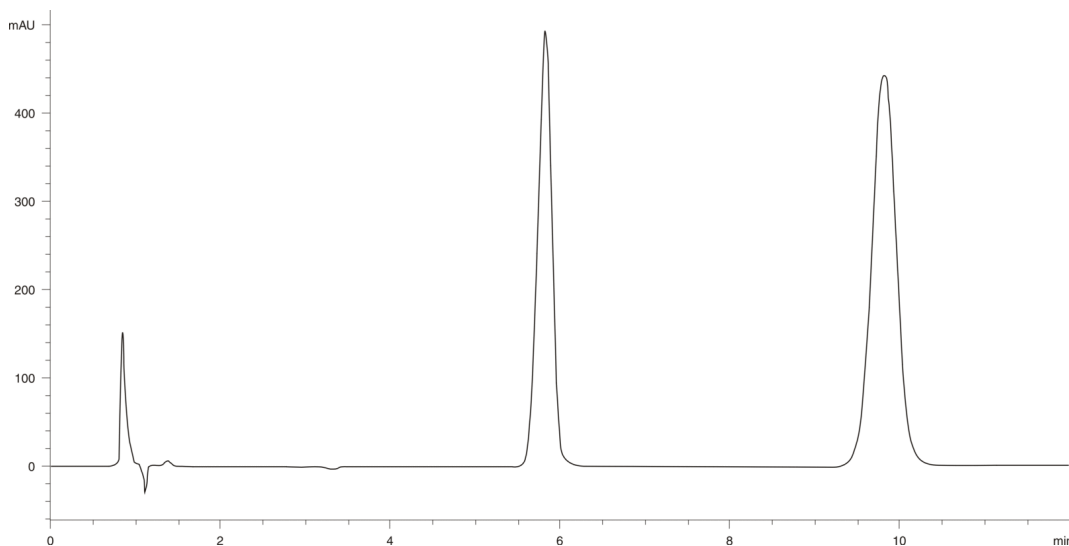

ANALYTES:

- 1 Felbamate
- 2 internal standard

