

Cough and Cold

Column: Prodigy™ 5 µm ODS-3 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4097-E0

Elution Type: Isocratic

Eluent A: Methanol

Eluent B: 30 mM KH₂PO₄, pH 3.0 with H₃PO₄

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	65	35	

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy™

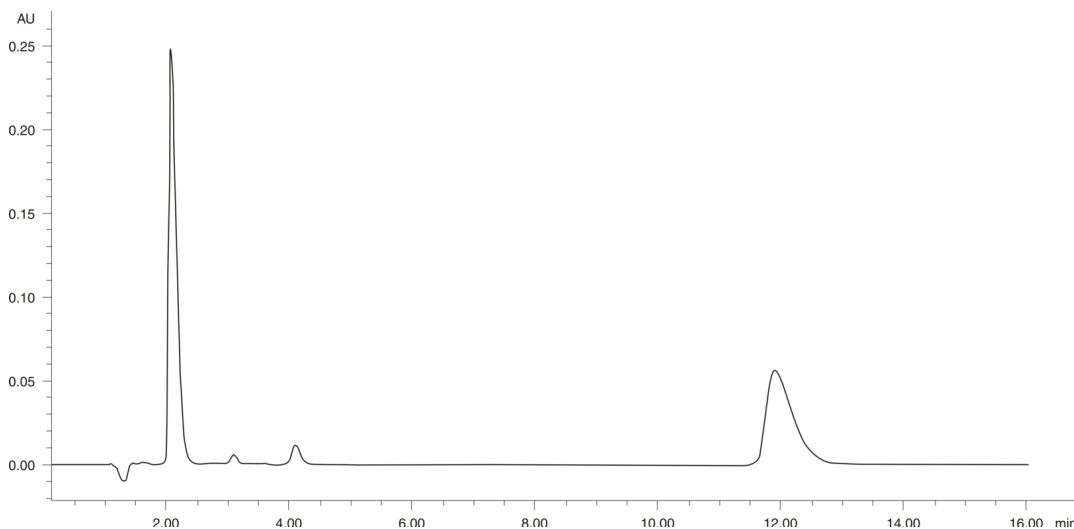

ANALYTES:

- 1 Pseudoephedrine
- 2 Methscopolamine
- 3 Chlorpheniramine maleate

