

Bretylum Tosylate - Assay for Related Compounds

Column: Prodigy™ 5 µm Phenyl-3 (PH-3) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4298-E0

Elution Type: Isocratic

Eluent A: 0.01 M Sodium 1-cotanesulfonate

Eluent B: Acetonitrile

Eluent C: Glacial acetic acid

Eluent D: Triethylamine

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C	Pct D
	1	0	79	18.5	2	0.5

Flow Rate: 1.9 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 265 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351

- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy™

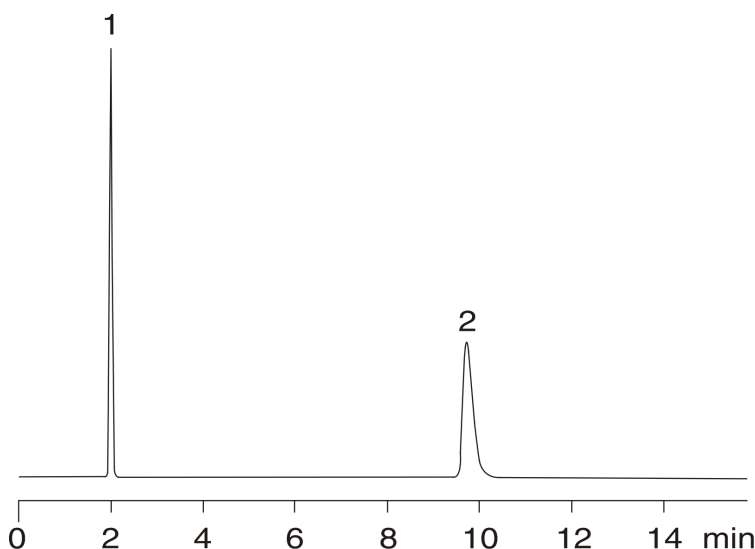

ANALYTES:

- 1 Tosylate ion
- 2 Bretylum

