

Food Additives

Column: Prodigy™ 5 µm Phenyl-3 (PH-3) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4298-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: 20 mM Potassium phosphate pH 3.0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	25	75

Flow Rate: 1.5 mL/min

Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)


Detector Info: Variable wavelength UV detector was changed during the run as follows:
0 to 1.6 min UV @ 230 nm
1.6 to 2 min UV @ 210 nm
>2 min UV @ 230 nm

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.
- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351
- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy™

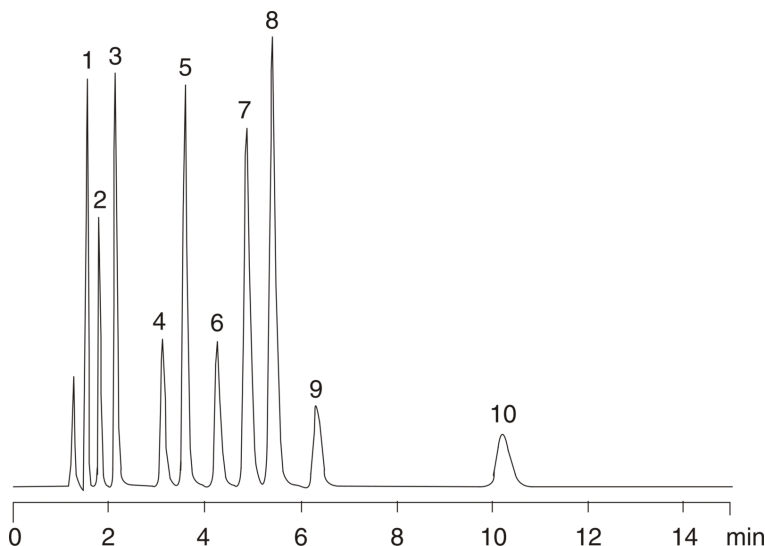

Food Additives

ANALYTES:

- 1 Saccharin
- 2 Aspartame
- 3 p-Hydroxybenzoic acid
- 4 Caffeine
- 5 Benzoic acid
- 6 p-Hydroxybenzoic acid methyl ester
- 7 p-Toluic acid
- 8 Dehydroacetic acid
- 9 p-Hydroxybenzoic acid ethyl ester
- 10 p-Hydroxybenzoic acid propyl ester

