

Benzoic Acids

Column: Prodigy™ 5 µm Phenyl-3 (PH-3) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4298-E0

Elution Type: Isocratic

Eluent A: 20 mM Potassium phosphate pH 2.5

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351

- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy™

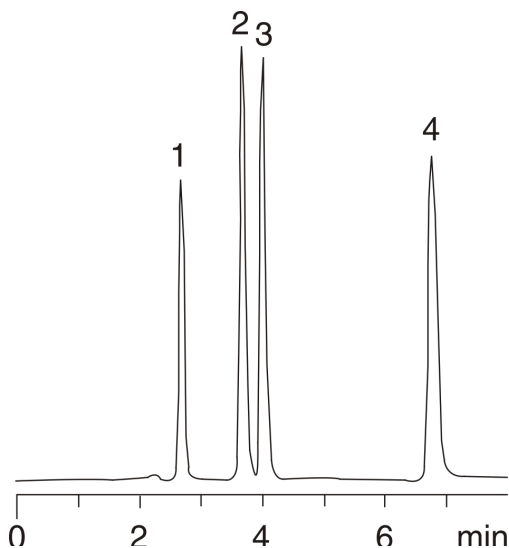

ANALYTES:

- 1 3,4-Dihydroxybenzoic acid
- 2 2,5-Dihydroxybenzoic acid
- 3 2,3-Dihydroxybenzoic acid
- 4 Benzoic acid

