

Coenzyme Q10

Column: Luna® 5 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: Ethanol

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	65	35	

Flow Rate: 1 mL/min


Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 275 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

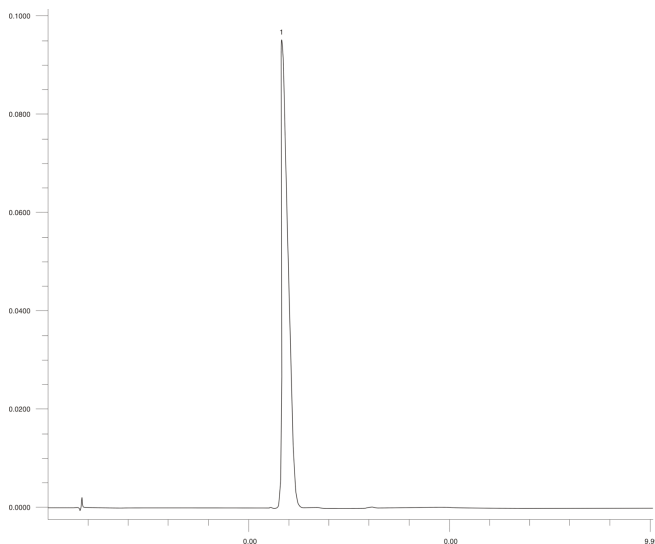

ANALYTES:

1 Coenzyme Q10



Products used in this application:

