


Test Conditions for Zebron ZB-1**Column:** Zebron[™] ZB-1, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G001-11**Oven Profile:** Isothermal @ 130°C**Carrier Gas:** Constant Flow Hydrogen, 1.93 mL/min**Injection:** Split 100:1 1 µL @ 250°C**Detection:** Flame Ionization (FID) (325°C)**Detector Info:** FID, 325°C, 20mL/min Nitrogen, 30mL/min**Products used in this application:****ANALYTES:**

- 1 Undecane
- 2 4-Chlorophenol
- 3 Tridecane
- 4 1-Undecanol
- 5 Dicyclohexylamine
- 6 Pentadecane

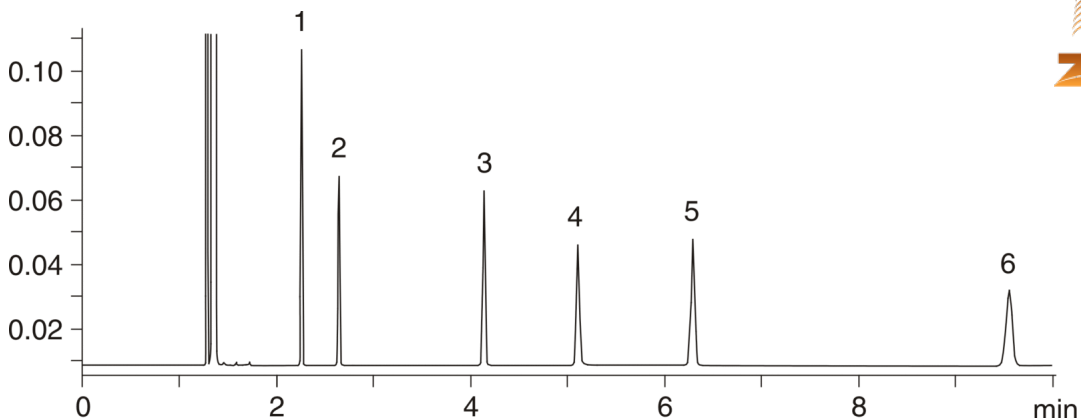