

St. John's Wort (*Hypericum perforatum*)

Column: Luna® 5 µm C18(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4252-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Methanol

Eluent C: 0.1N Ammonium Acetate, pH 9.7

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
1	0	50	30	20	

Flow Rate: 1.5 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 330 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

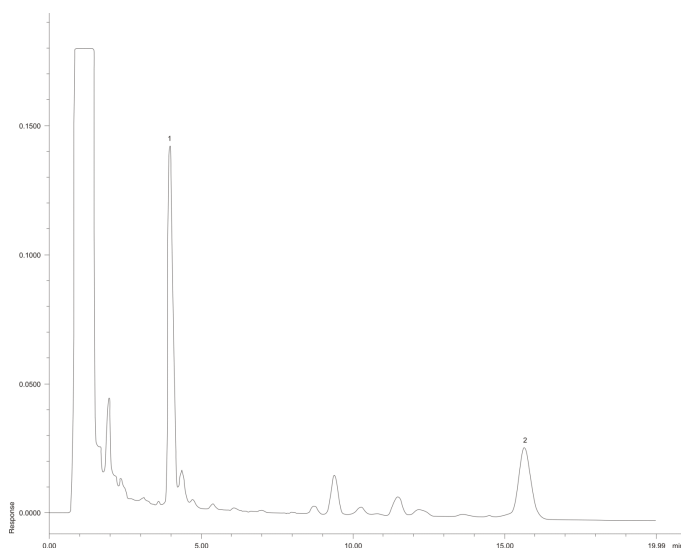

ANALYTES:

- 1 Pseudohypericin
- 2 Hypericin

