

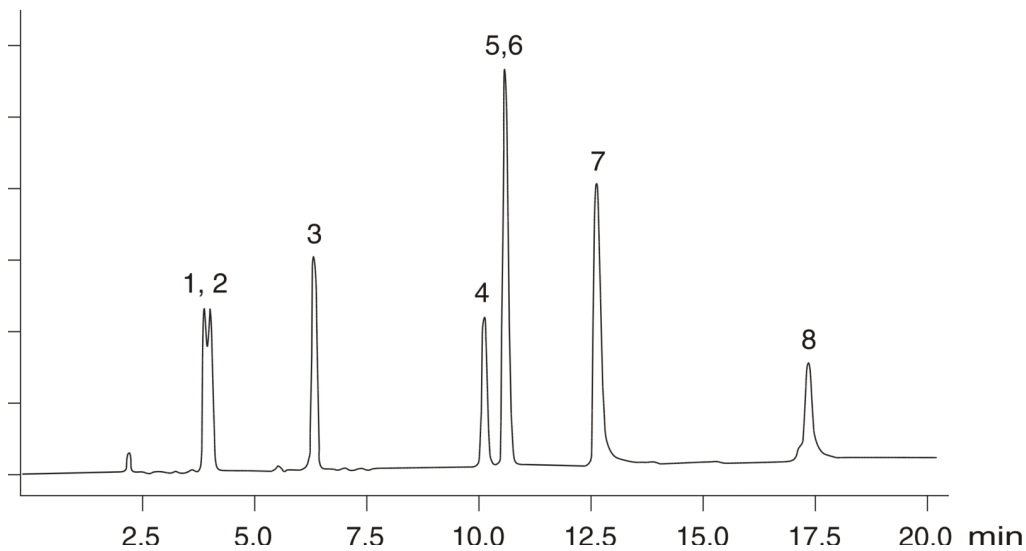
Antibacterials**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 20 mM KH₂PO₄, pH 2.5**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	5	75	25
	3	20	55	45

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Carbadox
- 2 Thiamphenicol
- 3 Furazolidone
- 4 Oxolinic acid
- 5 Sulfadimethoxine
- 6 Sulfaquinoxaline
- 7 Nalidixic acid
- 8 Piromidic acid

