

## Analytical-Like Performance

**Column:** Luna® 3 µm C18(2) 100 Å, LC Column 20 x 2 mm, Ea

**Dimensions:** 20 x 2 mm ID

**Order No:** 00M-4251-B0

**Elution Type:** Isocratic

**Eluent A:** Acetonitrile

**Eluent B:** Methanol

**Eluent C:** 20 mM Potassium phosphate, pH 7.0

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B | Pct C |
|-------------------|----------|------------|-------|-------|-------|
| 1                 | 0        | 40         | 30    | 30    |       |

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 2.0mm, 10/Pk Part No.: AJ0-4286

- Holder Part No.: KJ0-4282

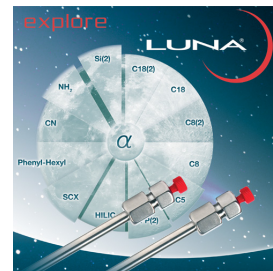

Products used in this application:

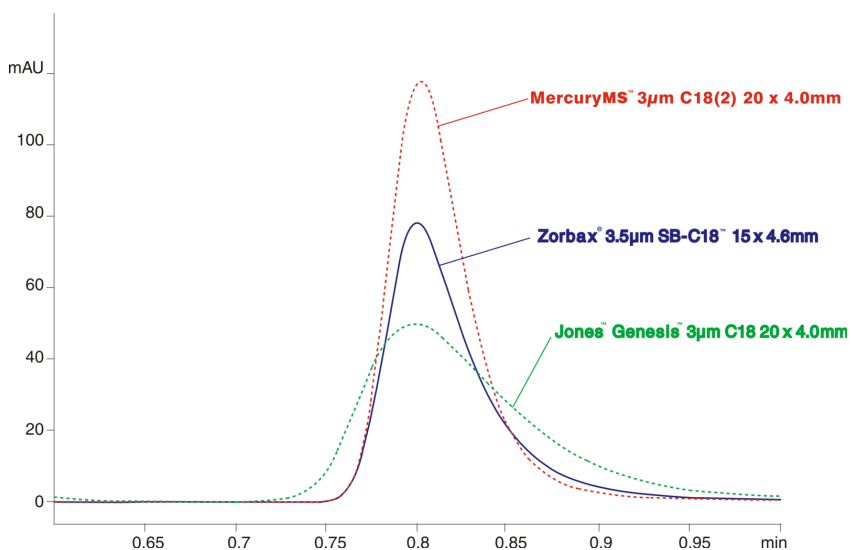

## ANALYTES:

1 Nortriptyline

