

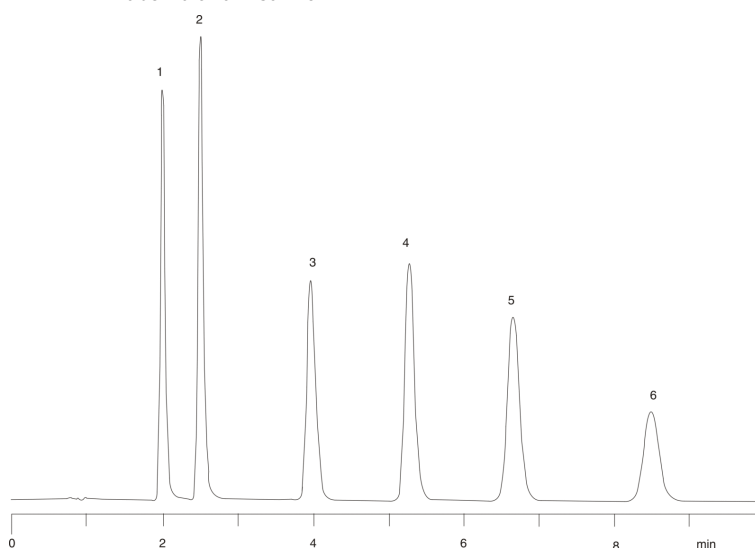
Tricyclic Antidepressants**Column:** SynergiTM 4 μ m Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4337-E0**Elution Type:** Isocratic**Eluent A:** 20 mM Potassium phosphate pH 7**Eluent B:** Methanol**Eluent C:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
	1	0	30	35	35

Flow Rate: 2 mL/min**Col. Temp.:** 30 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Desipramine
- 2 Nortriptyline
- 3 Doxepin
- 4 Imipramine
- 5 Amitriptyline
- 6 Clomipramine

