

Nicotine and Related Compounds

Column: Synergi™ 4 µm Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4337-E0

Elution Type: Gradient

Eluent A: 10 mM Ammonium acetate in water

Eluent B: 10 mM Ammonium acetate in methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	0.25	90	10
	3	3.75	0	100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

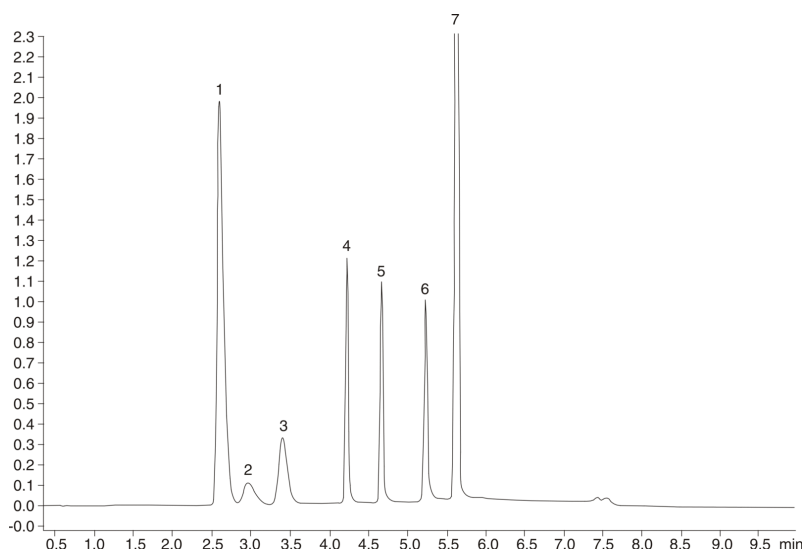

ANALYTES:

- 1 Cotinine-N'-oxide
- 2 1'(R)-2'(S)-cis-Nicotine-N'-oxide
- 3 1'(S)-2'(S)-trans-Nicotine-N'-oxide
- 4 Hydroxycotinine
- 5 Cotinine
- 6 Nicotine
- 7 Picric acid

