

Metal Sensitivity Test - Luna 5 μ C5

Column: Luna[®] 5 μ m C5 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4043-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	50	50	

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)


Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C5 4 x 3.0mm, 10/Pk Part No.: AJ0-4293

- Holder Part No.: KJ0-4282



Products used in this application:

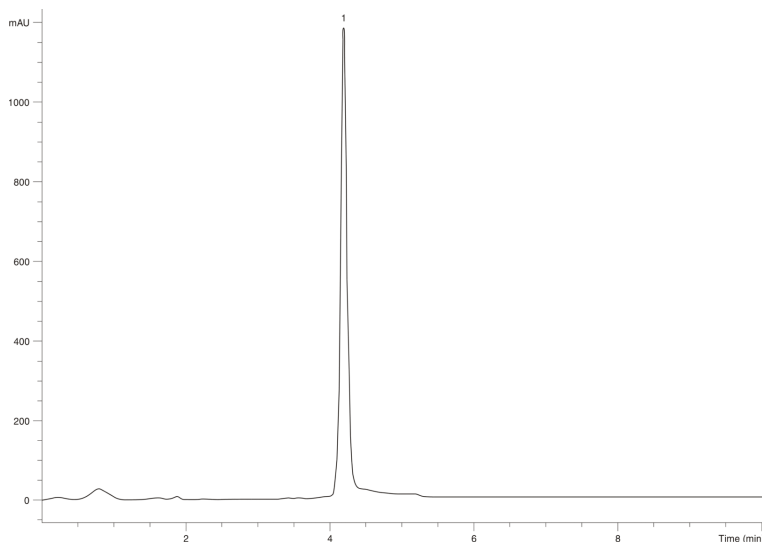

ANALYTES:

1 8-Hydroxyquinoline

