

Increase Sample Through-put 20x4.0mm

Column: Luna® 3 µm C18(2) 100 Å, LC Column 20 x 4 mm, Ea

Dimensions: 20 x 4 mm ID

Order No: 00M-4251-D0

Elution Type: Gradient

Eluent A: Water w/0.1% formic acid

Eluent B: Acetonitrile w/0.1% formic acid

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------

Profile: **1** **0** **95** **5**

2 2.4 5 95

Flow Rate: 1.6 mL/min

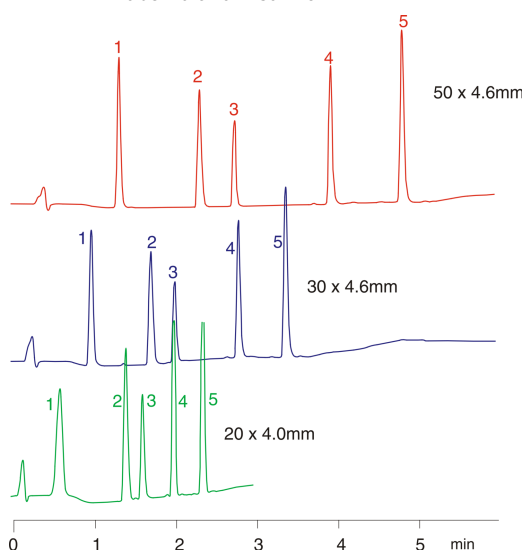
Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Acetaminophen
- 2 Propranolol
- 3 Imipramine
- 4 Naproxen
- 5 Valerophenone



Products used in this application:

