

pH Flexibility - Figure 7d

Column: Luna® 5 µm C18(2) 100 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4252-E0

Elution Type: Gradient

Eluent A: Water w/10mM ammonium acetate, pH 5.0

Eluent B: Acetonitrile w/10mM ammonium acetate, pH 5.0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	6	5	95

Flow Rate: 2 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

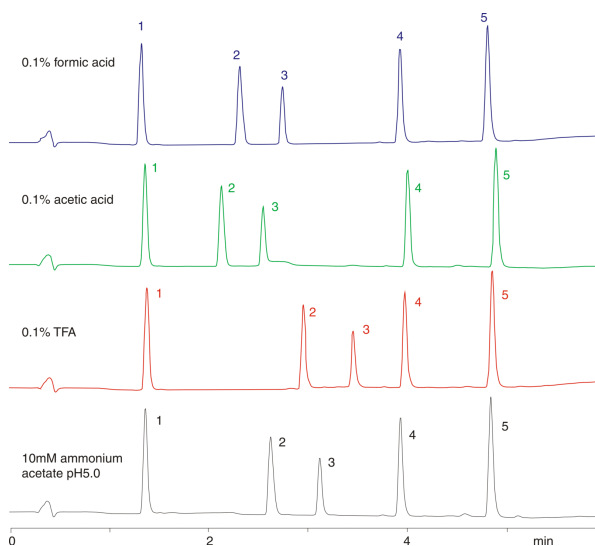

ANALYTES:

- 1 Acetaminophen
- 2 Propranolol
- 3 Imipramine
- 4 Naproxen
- 5 Valerophenone

