

CNS Stimulant

Column: Luna® 5 µm Phenyl-Hexyl 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4257-E0

Elution Type: Isocratic

Eluent A: Methanol

Eluent B: 0.15M Ammonium acetate, pH 4.0

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	6	94

Flow Rate: 1.5 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 257 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351

- Holder Part No.: KJ0-4282



Products used in this application:

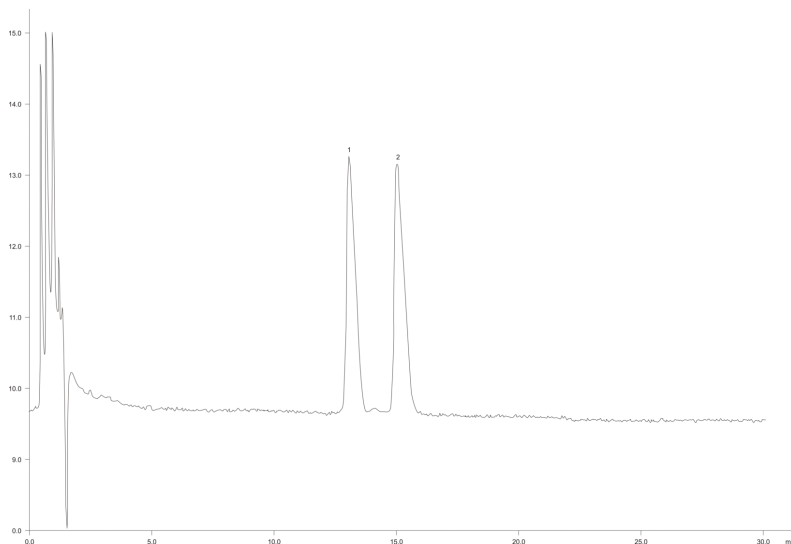

ANALYTES:

- 1 Ephedrine
- 2 Pseudoephedrine

