


**Ginkgo Leaf (Ginkgo biloba)****Column:** Zebtron<sup>™</sup> ZB-1, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G001-11**Oven Profile:** 120°C to 340°C @ 5°C/min**Carrier Gas:** Constant Pressure Helium, 1.7 mL/min**Injection:** Splitless :1 1 µL @ 280°C**Detection:** Flame Ionization (FID) (300°C)

1122

**Products used in this application:****ANALYTES:**

- 1 Benzophenone
- 2 Bilobalide
- 3 Ginkgolide A
- 4 Ginkgolide B

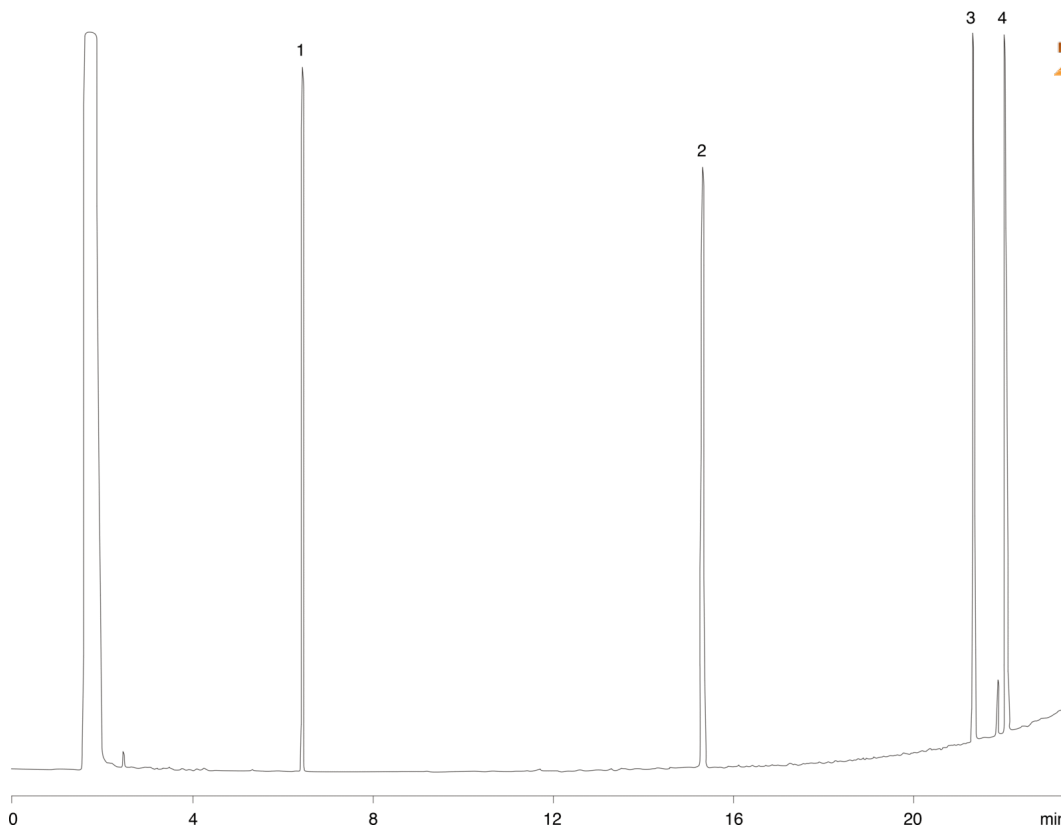