

## Chinese Green Tea (Comellia sinensis)

**Column:** Prodigy™ 5 µm ODS-3 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4097-E0

**Elution Type:** Isocratic

**Eluent A:** 0.2% Phosphoric Acid

**Eluent B:** Acetonitrile

**Eluent C:** Tetrahydrofuran

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B | Pct C |
|-------------------|----------|------------|-------|-------|-------|
| 1                 | 1        | 0          | 86    | 12.5  | 1.5   |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 25 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

**Prodigy™**

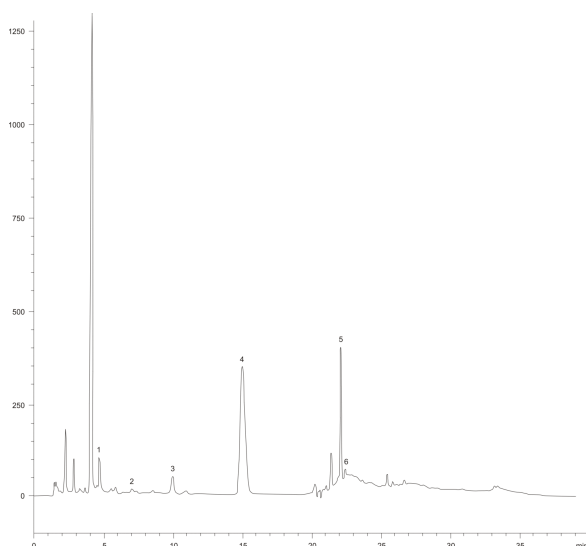

### ANALYTES:

- 1 Epigallocatechin
- 2 Catechin
- 3 Epicatechin
- 4 Epigallocatechin gallate
- 5 Epicatechin gallate
- 6 Catechin gallate

