

Sulfa Drugs

Column: Luna® 5 µm C5 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4043-E0

Elution Type: isocratic

Eluent A: Acetonitrile

Eluent B: Water, pH 9.0, pH 9

Gradient	Step No.	Time (min)	Pct A	Pct B
----------	----------	------------	-------	-------

Profile: 1 0 3 97

Flow Rate: 1 mL/min

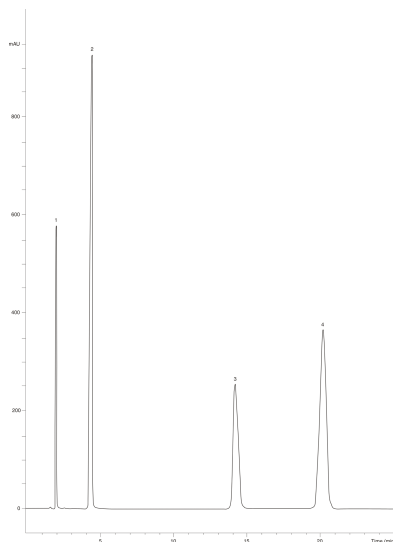
Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 270 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C5 4 x 3.0mm, 10/Pk Part No.: AJ0-4293

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1** Sulfanilic acid
- 2** Sulfanilamide
- 3** Sulfathiazole
- 4** Sulfapyridine



Products used in this application:

