

Quinones

Column: Luna® 3 µm C8(2) 100 Å, LC Column 75 x 4.6 mm, Ea

Dimensions: 75 x 4.6 mm ID

Order No: 00C-4248-E0

Elution Type: Isocratic

Eluent A: Methanol

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	80	20	

Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 251 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282



Products used in this application:

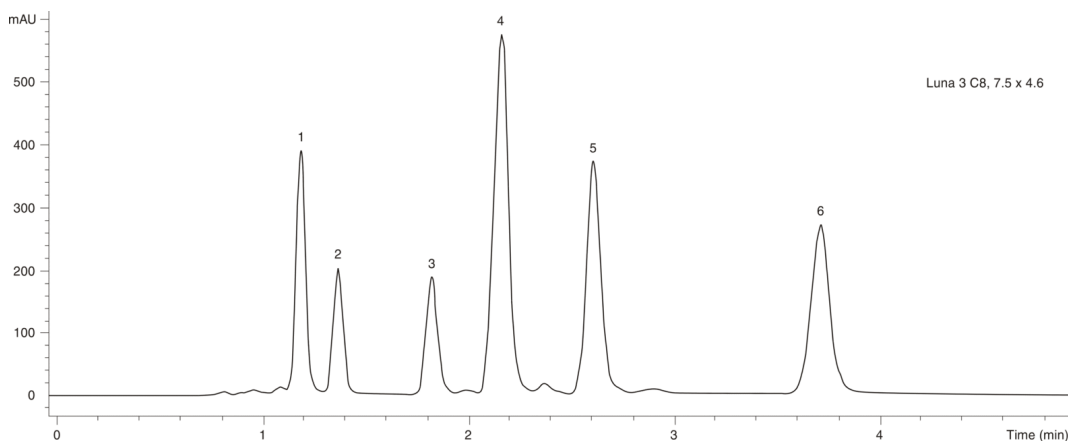

ANALYTES:

- 1 1,4-Naphthoquinone
- 2 2-Methyl-1,4-naphthquinone
- 3 9,10-Anthraquinone
- 4 1,5-Dichloro-9,10-anthraquinone
- 5 2-Methyl

