

Food Additives on Luna 5u Phenyl-Hexyl

Column: Luna® 5 µm Phenyl-Hexyl 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4257-E0

Elution Type: Gradient

Eluent A: 50 mM KH₂PO₄ + 0.1% H₃PO₄

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	18	25	75
	3	30	25	75

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-4351

- Holder Part No.: KJ0-4282



Products used in this application:

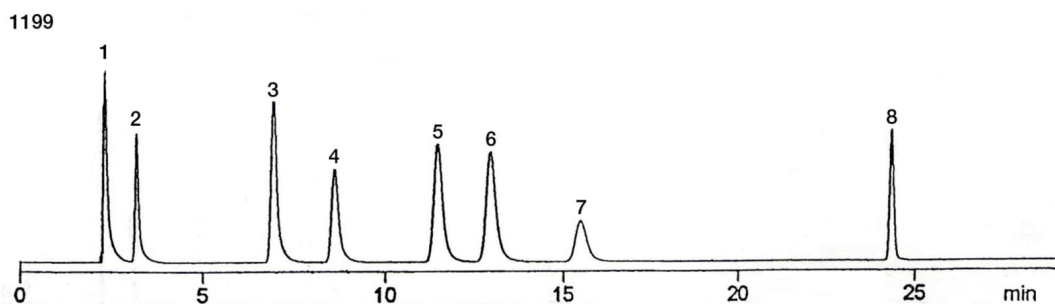

ANALYTES:

- 1 Saccharin
- 2 p-Hydroxybenzoic acid
- 3 Sorbic acid
- 4 p-Hydroxybenzoic acid methyl ester
- 5 Dehydroacetic acid
- 6 p-Toluic acid
- 7 p-Hydroxybenzoic acid ethyl ester
- 8 n-Propyl p-hydroxybenzoate

