

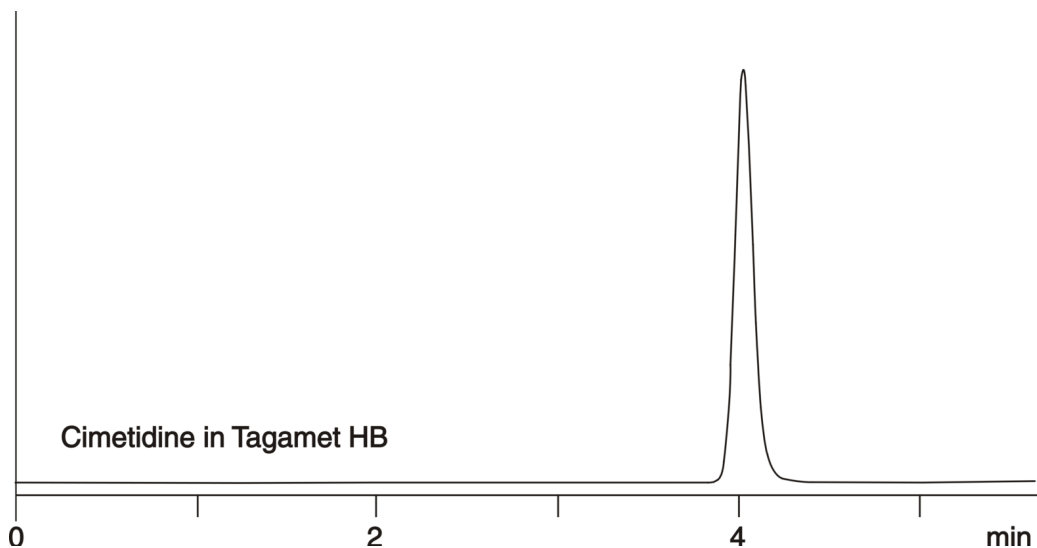
Cimetidine**Column:** SynergiTM 4 μ m Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4337-E0**Elution Type:** Isocratic**Eluent A:** 20mM Potassium phosphate pH 7**Eluent B:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	60	40

Flow Rate: 1 mL/min**Col. Temp.:** 30 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**ANALYTES:**

1 Cimetidine

