


Gasoline Range Organics**Column:** ZebtronTM ZB-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G001-11**Oven Profile:** 35°C (2 min) to 110°C at 6°C/min**Carrier Gas:** Constant Pressure Helium, 0.88 mL/min**Injection:** Split 20:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Sample analytes 200ppm GROs in methanol**Products used in this application:****ANALYTES:**

- 1 3-Methylpentane
- 2 2,2,4-Trimethylpentane
- 3 Benzene
- 4 Toluene
- 5 Ethylbenzene
- 6 m-Xylene
- 7 o-Xylene
- 8 1,2,4-Trimethylbenzene
- 9 Naphthalene

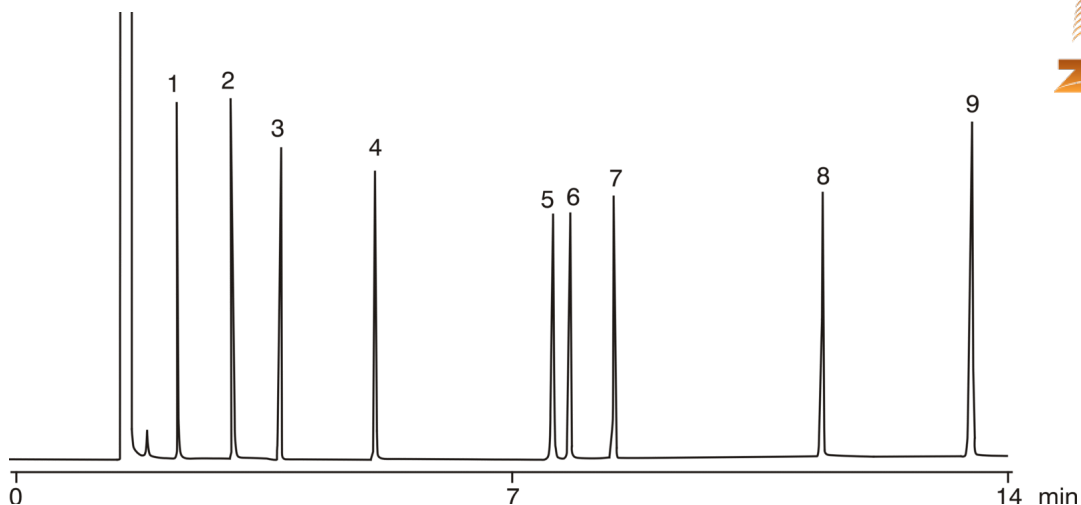