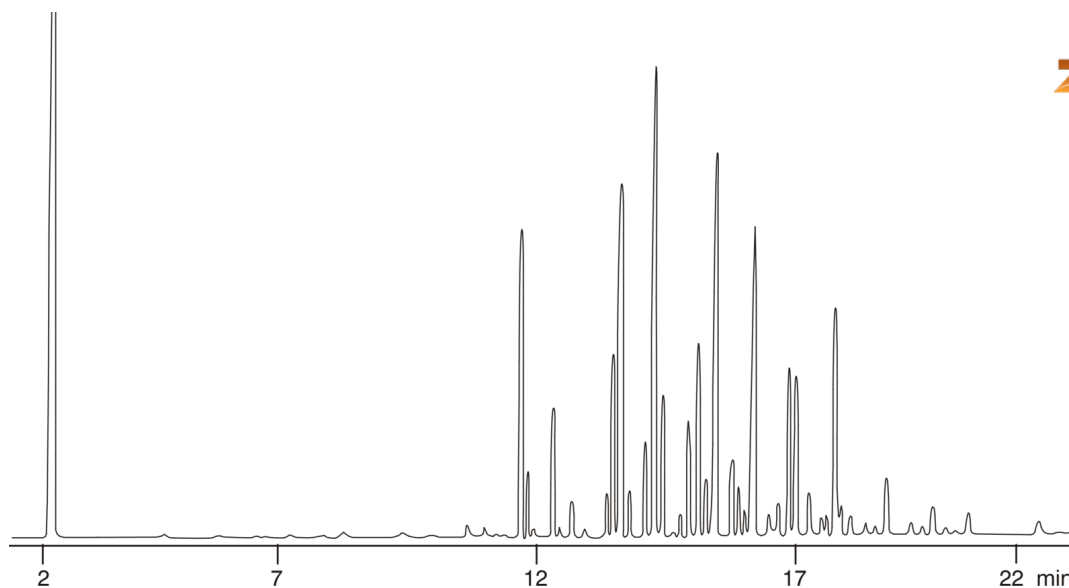


Aroclor 1254: EPA Method 8082**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G002-11**Oven Profile:** 160°C to 260°C at 6°C/min (hold 4 min)**Carrier Gas:** Constant Flow Helium, 0.88 mL/min**Injection:** Split 100:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Sample analytes 100ppm aroclor 1254 in hexane**Products used in this application:****ANALYTES:****1** Arochlor 1254