

Aldehydes

Column: Luna® 5 µm C18(2) 100 Å, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4252-B0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water, 0.1% H₃PO₄

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	50	50	

Flow Rate: 0.35 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 370 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 2.0mm, 10/Pk Part No.: AJ0-4286

- Holder Part No.: KJ0-4282



Products used in this application:

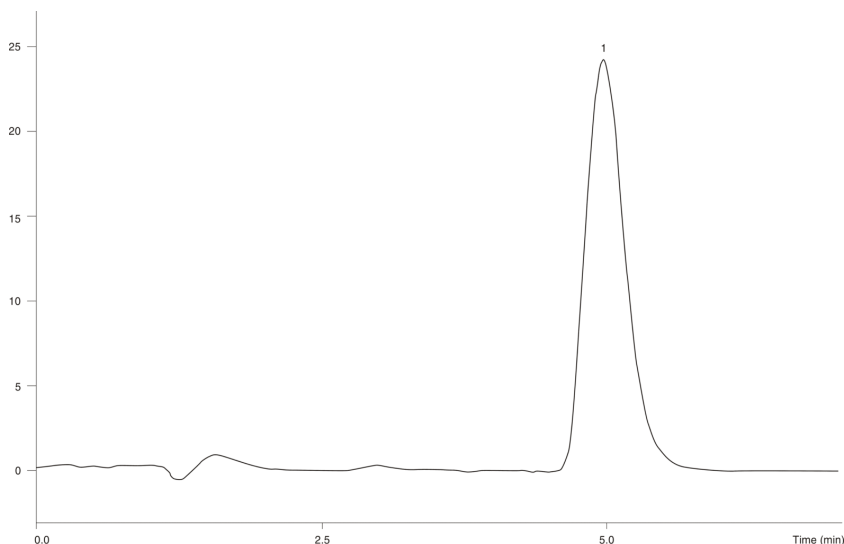

ANALYTES:

1 DNPH formaldehyde

