

Flavors

Column: Luna® 3 µm Phenyl-Hexyl 100 Å, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4256-B0

Elution Type: Isocratic

Eluent A: 0.15% TFA in Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	63	37	

Flow Rate: 0.7 mL/min

Col. Temp.: 60 °C

Detection: Mass Spectrometer @ (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Phenyl 4 x 2.0mm, 10/Pk Part No.: AJ0-4350

- Holder Part No.: KJ0-4282



Products used in this application:

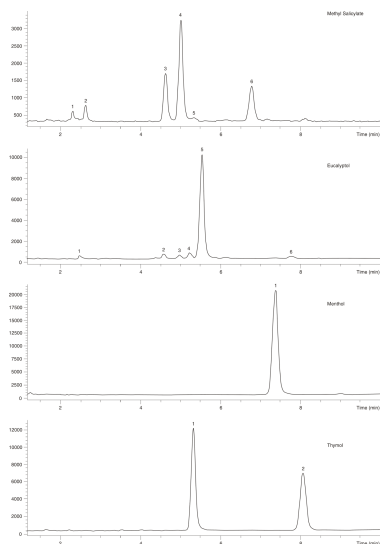

ANALYTES:

- 1 Methyl salicylate
- 2 Eucalyptol
- 3 Menthol
- 4 Thymol

