

Carbamate/Urea Pesticides

Column: Luna® 3 µm CN 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4254-E0

Elution Type: Gradient

Eluent A: Hexane

Eluent B: 80:20 Methylene chloride:Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	2	85	15

Flow Rate: 5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282

12945



Products used in this application:

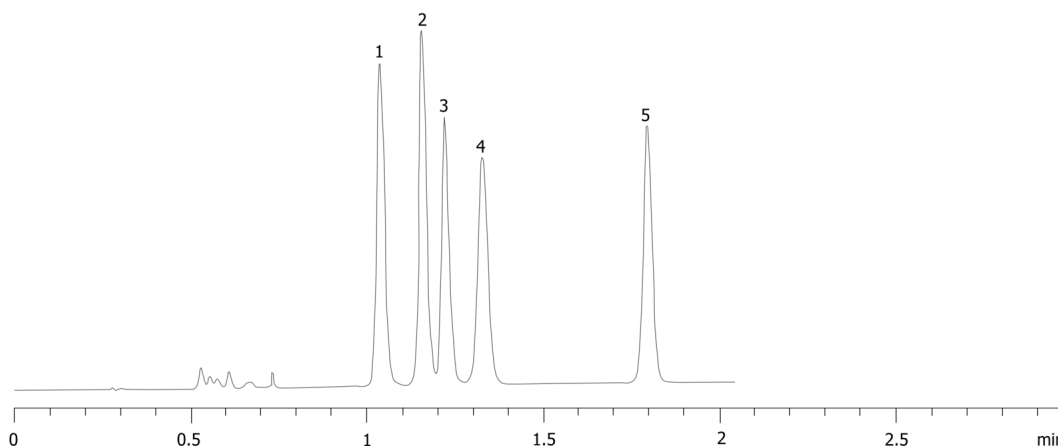

ANALYTES:

- 1 Carbofuran
- 2 Barban
- 3 Carbaryl
- 4 Tebuthiuron
- 5 Diuron

