

Plasticizers

Column: Luna® 3 µm CN 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4254-E0

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: 80:20 Methylene chloride:Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	85	15	

Flow Rate: 4 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282

12948



Products used in this application:

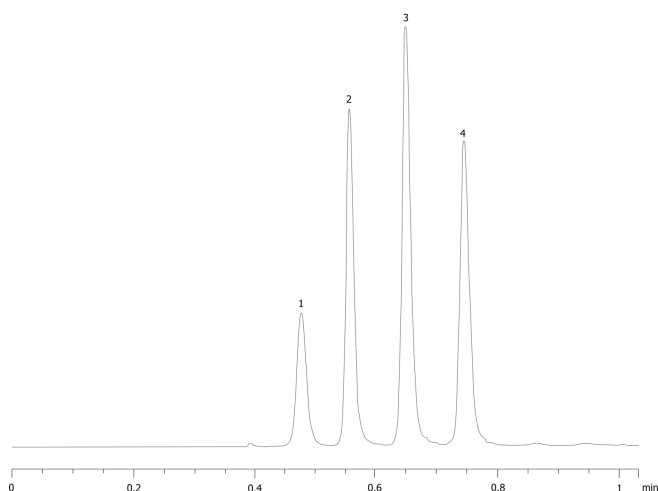

ANALYTES:

- 1 Diisooctyl phthalate
- 2 Dibutyl phthalate
- 3 Diethyl phthalate
- 4 Dimethyl phthalate

