

Cyano Comparisons - Carbamate/Urea Pesticides and Herbicides

Column: Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4255-E0

Elution Type: Gradient

Eluent A: Acetonitrile

Eluent B: Water

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	18	82
	2	5	65	35

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282

12951



Products used in this application:

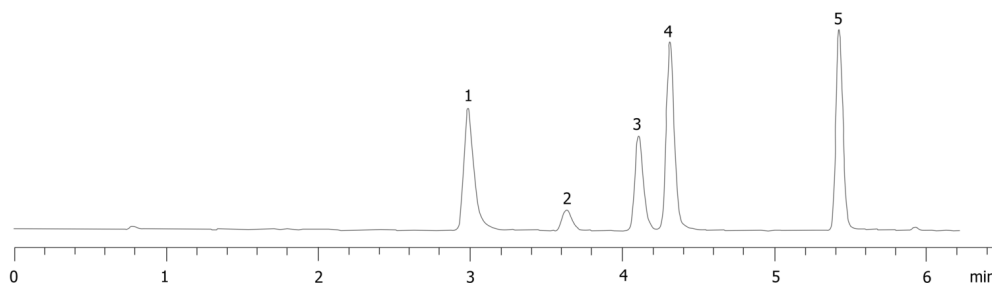

ANALYTES:

- 1 Tebuthiuron
- 2 Carbofuran
- 3 Carbaryl
- 4 Diuron
- 5 Barban

