

## Cyano Comparisons - Rhodamine Dyes

**Column:** Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4255-E0

**Elution Type:** Isocratic

**Eluent A:** 80:20 20 mM Potassium phosphate pH 3.0:Methanol

**Eluent B:** Acetonitrile

**Gradient**

| Step No. | Time (min) | Pct A | Pct B |
|----------|------------|-------|-------|
| 1        | 0          | 50    | 50    |

**Profile:**

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282

12953



Products used in this application:

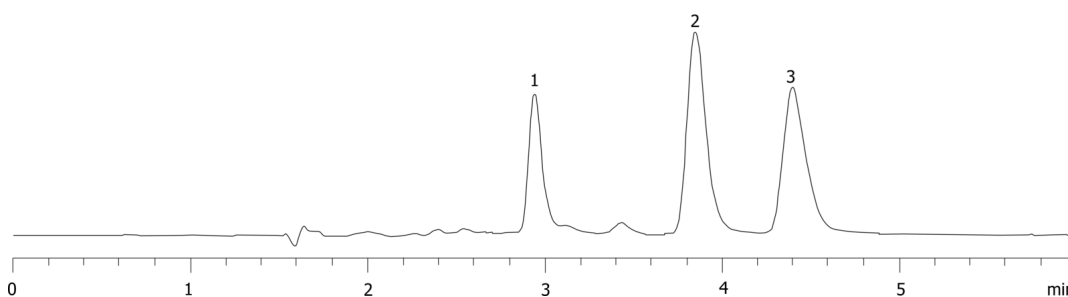

### ANALYTES:

- 1 Rhodamine 123
- 2 Rhodamine B
- 3 Rhodamine 6G

