


Herbicides: EPA Method 8051**Column:** Zebtron[™] ZB-1701, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 14% Cyanopropylphenyl 86% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G006-11**Oven Profile:** 130°C to 280°C at 12°C/min (hold 10 min)**Carrier Gas:** Constant Flow Helium, 0.88 mL/min**Injection:** Split 10:1 1 μ L @ 250°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** 100ppm herbicide mixture in ethyl acetate**Products used in this application:**

Herbicides: EPA Method 8051

ANALYTES:

- 1 Eptam[®]
- 2 Sutan[®]
- 3 Tillam[®]
- 4 Molinate
- 5 Cycloate
- 6 Treflan[®]
- 7 Atrazine
- 8 Terbacil
- 9 Sencor[®]
- 10 Bromacil
- 11 Isopropalin
- 12 Oxyfluorfen
- 13 Velpar[®]

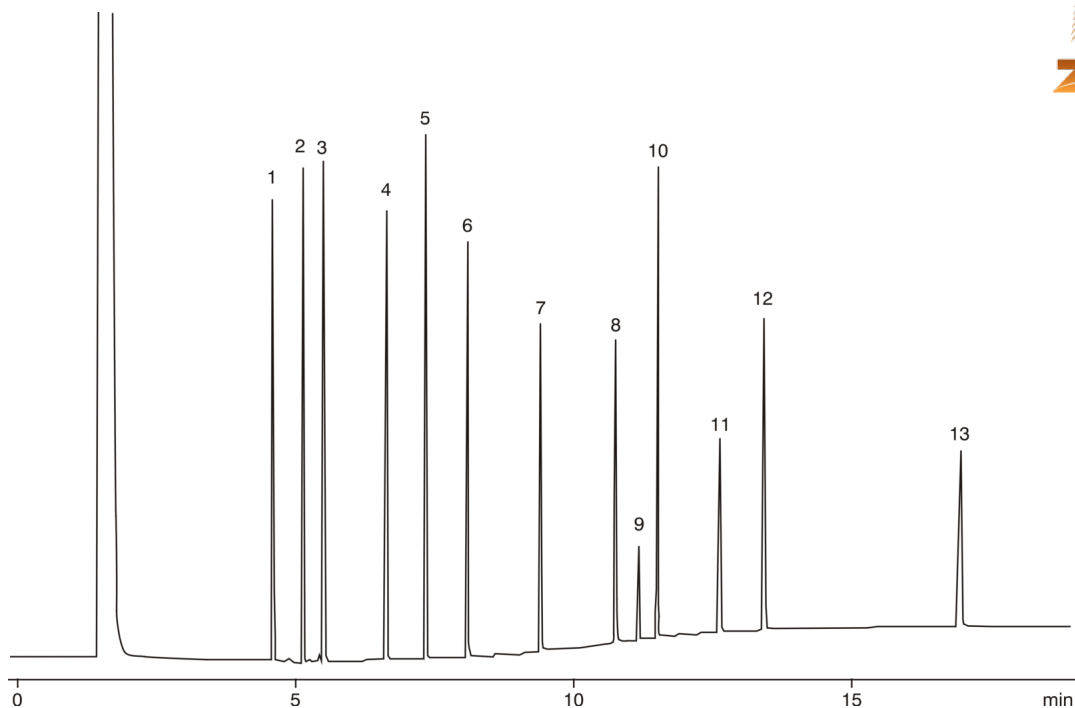