


Fragrance Analysis

Column: Zebron™ ZB-WAX, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: Polyethylene Glycol (PEG)
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G007-11
Oven Profile: 50°C (8 min) to 250°C @ 25°C/min (hold 2 min)

Carrier Gas: Constant Flow Helium, 0.88 mL/min
Injection: Split 50:1 1 µL @ 250°C
Detection: Flame Ionization (FID) (300°C)
Analyst Note: Sample 0.5% Early Spring fragrance in methanol



Products used in this application:

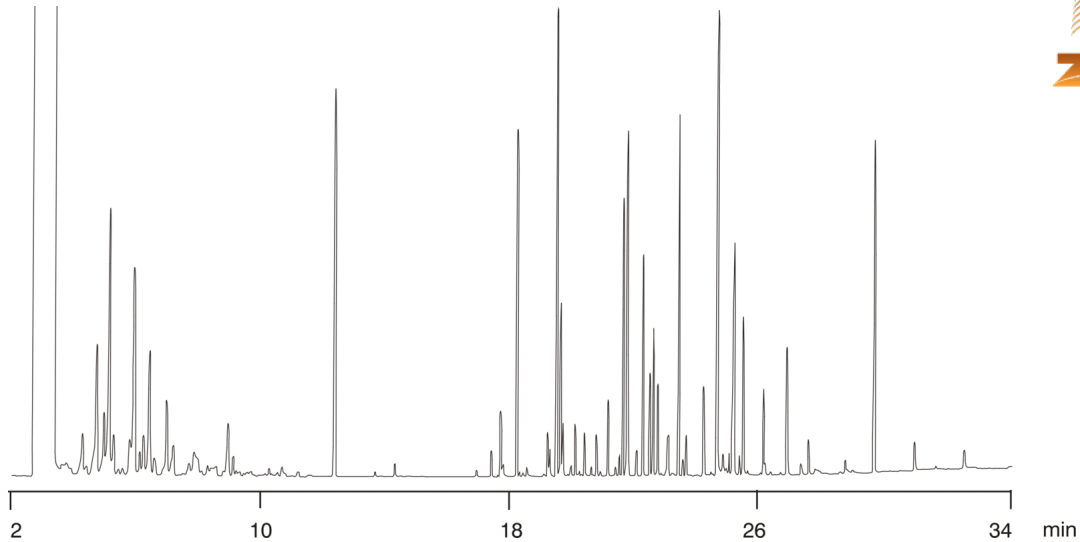

ANALYTES:

1 unidentified

