

Organic Compounds

Column: Luna® 3 µm C18(2) 100 Å, LC Column 75 x 4.6 mm, Ea

Dimensions: 75 x 4.6 mm ID

Order No: 00C-4251-E0

Elution Type: Isocratic

Eluent A: 40% Acetonitrile + 0.1% phosphoric acid

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 240 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

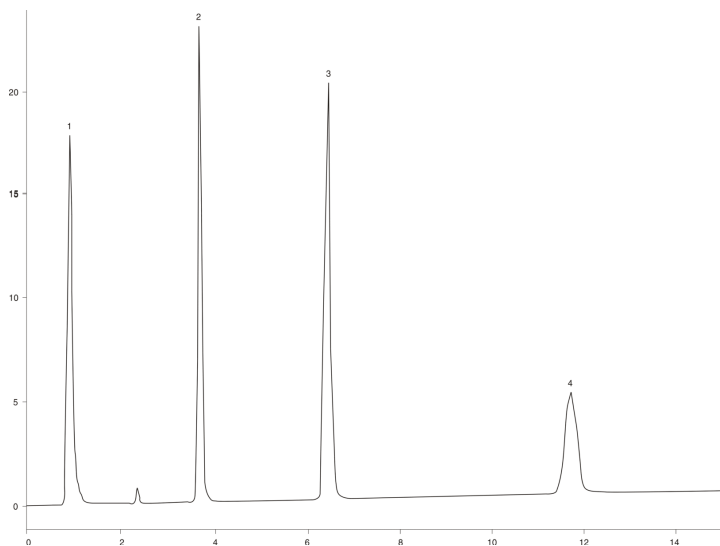

ANALYTES:

- 1 Ammonium nitrate
- 2 Ethyl paraben
- 3 Butyl paraben
- 4 Butyl-p-amino ben

