

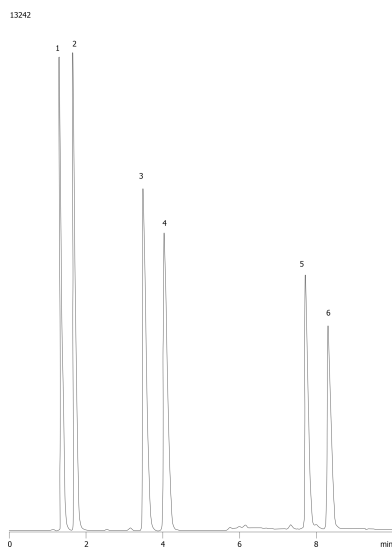
Environmental**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Water**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	30	70
	2	4	30	70
	3	10	0	100

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1** Dimethylphthalate
- 2** Diethylphthalate
- 3** Di-n-butyl phthalate
- 4** Butylbenzyl phthalate
- 5** bis(2-Ethylhexyl)phthalate
- 6** Di-n-octylphthalate

