

Antibiotic

Column: Prodigy™ 5 µm ODS-3 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4097-E0

Elution Type: Isocratic

Eluent A: Acetonitrile/0.05 M NaH₂PO₄ (20:80)

Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:

Flow Rate: 0.25 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

Prodigy™

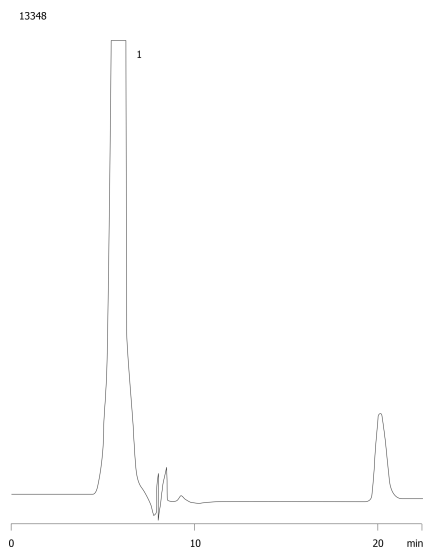

ANALYTES:

1 Ketamine

