

Acebutolol

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Hexane/Dichloroethane/Ethanol-TFA (55:35:10), with ethanol-TFA (20:1)

Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 244 nm (ambient)

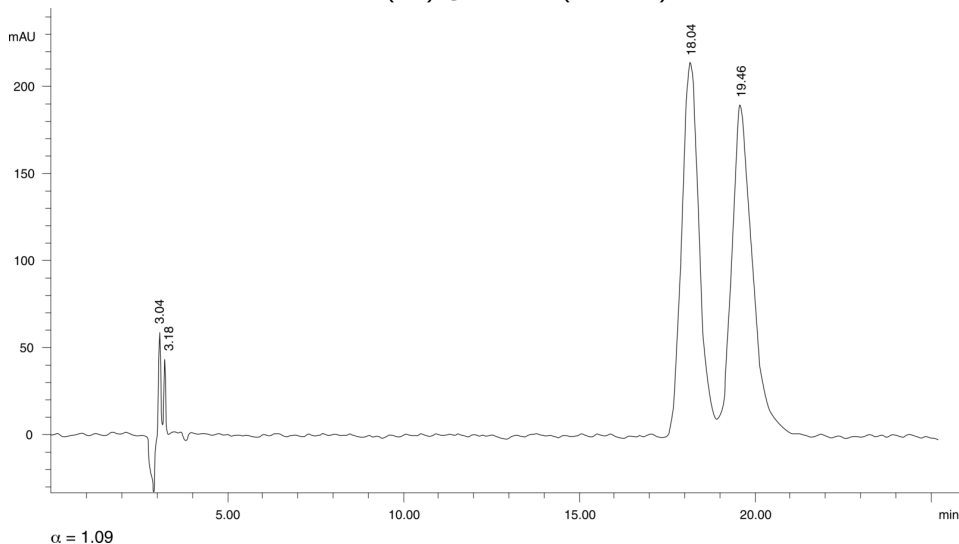

ANALYTES:

1 Acebutolol

Retention Time: 18.04 min

2 Acebutolol

Retention Time: 19.46 min

