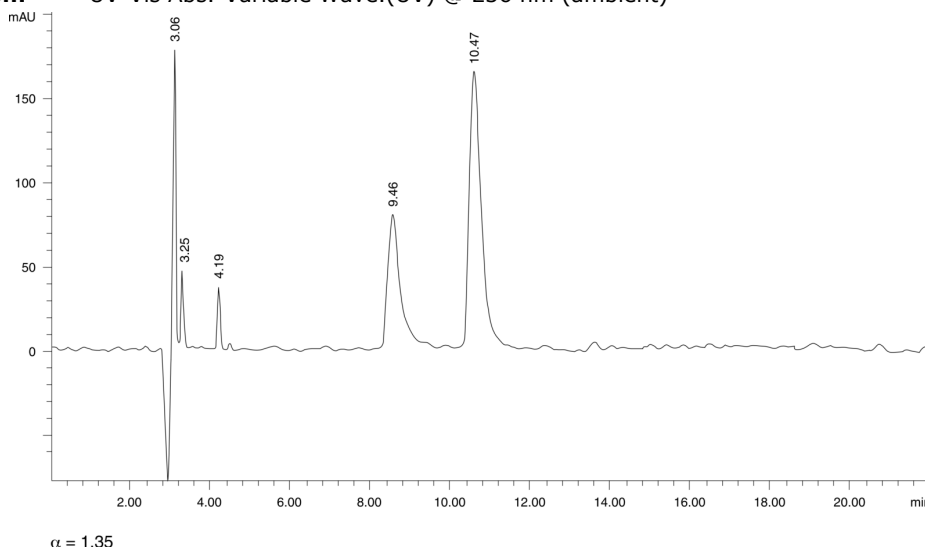


Clemastine**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (58:35:7) w/ethanol-TFA (20:1)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 236 nm (ambient)**ANALYTES:****1** Clemastine

Retention Time: 8.46 min

2 Clemastine

Retention Time: 10.47 min

