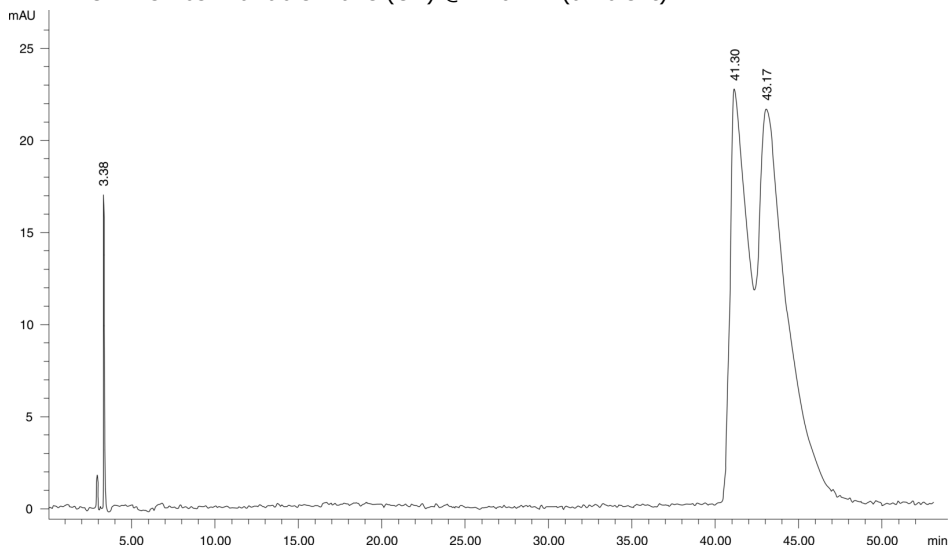


Isoxsuprine

Column: N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (60:35:5) ethanol-TFA 20:1

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 276 nm (ambient)

ANALYTES:

1 Isoxsuprine

Retention Time: 41.3 min

2 Isoxsuprine

Retention Time: 43.17 min

