


Labetalol


Column: Chirex® (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3020-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (60:35:5) ethanol-TFA 20:1
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 308 nm (ambient)

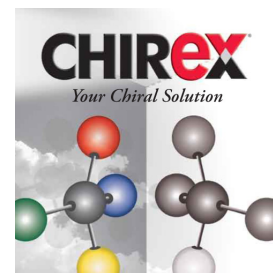

Products used in this application:

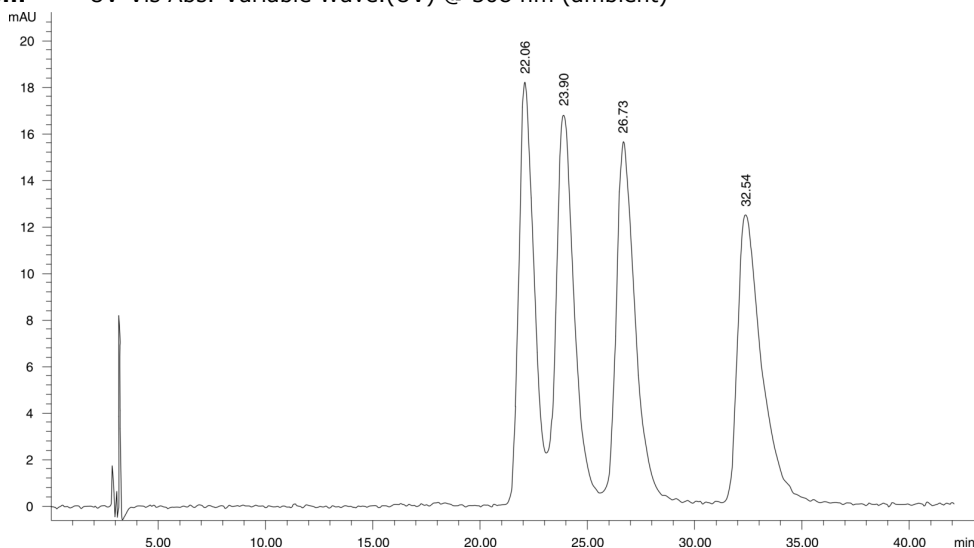CHIREX®

ANALYTES:

- 1 Labetalol
Retention Time: 22.06 min
- 2 Labetalol
Retention Time: 23.9 min
- 3 Labetalol
Retention Time: 26.73 min
- 4 Labetalol
Retention Time: 32.54 min

