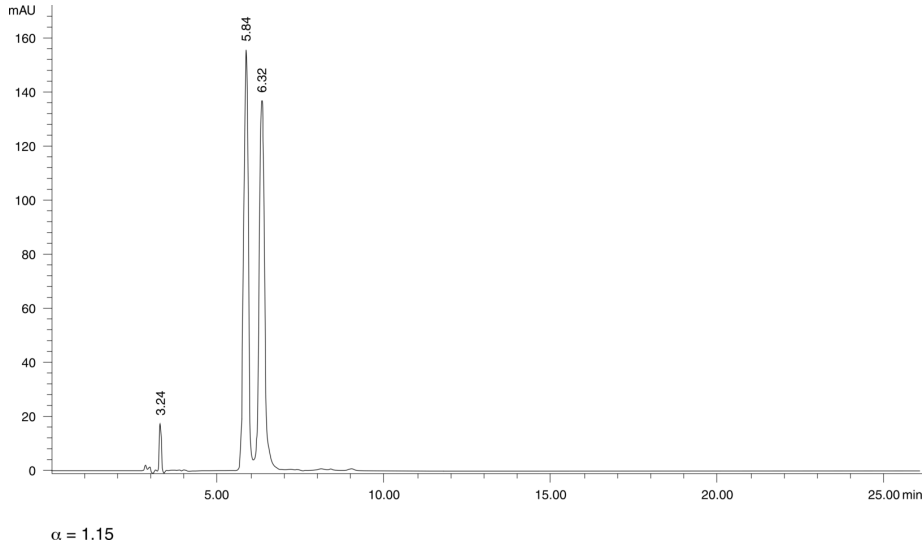


Methoxyverapamil

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (55:35:10)
Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 278 nm (ambient)

**ANALYTES:****1** Methoxyverapamil

Retention Time: 5.84 min

2 Methoxyverapamil

Retention Time: 6.32 min

