

Midodrine

Column: Chirex® (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-3020-E0

Elution Type: Isocratic


Eluent A: Hexane/Dichloroethane/Ethanol-TFA (58:35:7)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 292 nm (ambient)

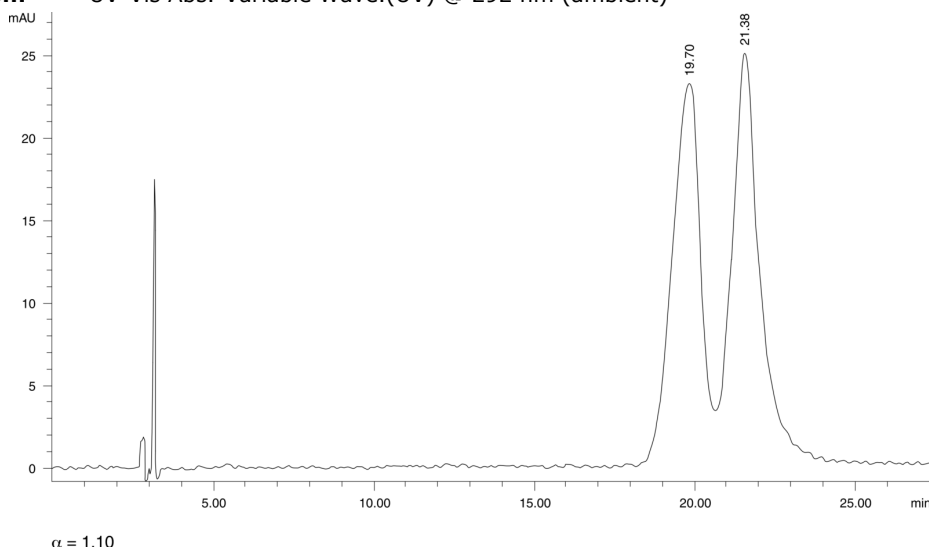

Products used in this application:

CHIREX®

ANALYTES:

1 Midodrine

Retention Time: 19.7 min

2 Midodrine

Retention Time: 21.38 min

