


Nicotine**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (55:35:10) ethanol-TFA 20:1**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)**ANALYTES:****1** Nicotine

Retention Time: 37.95 min

2 Nicotine

Retention Time: 39.83 min

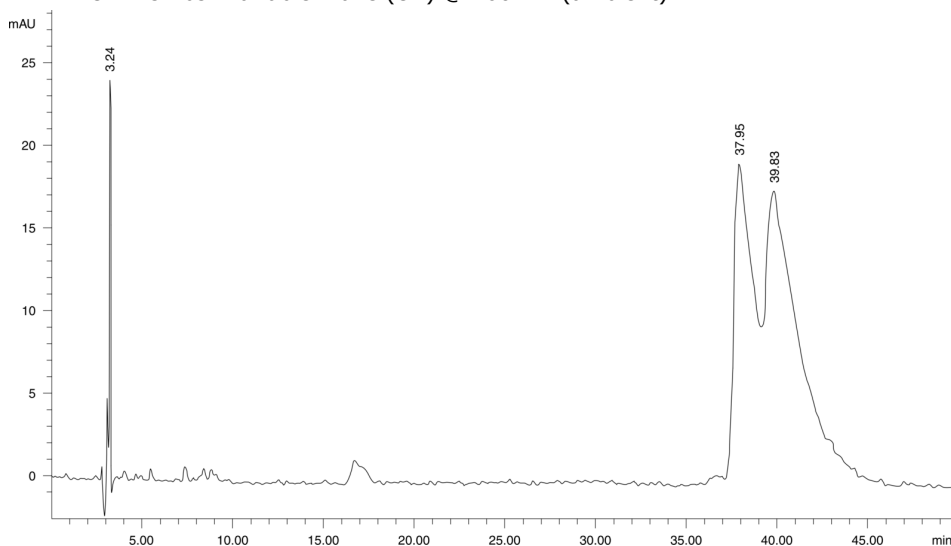