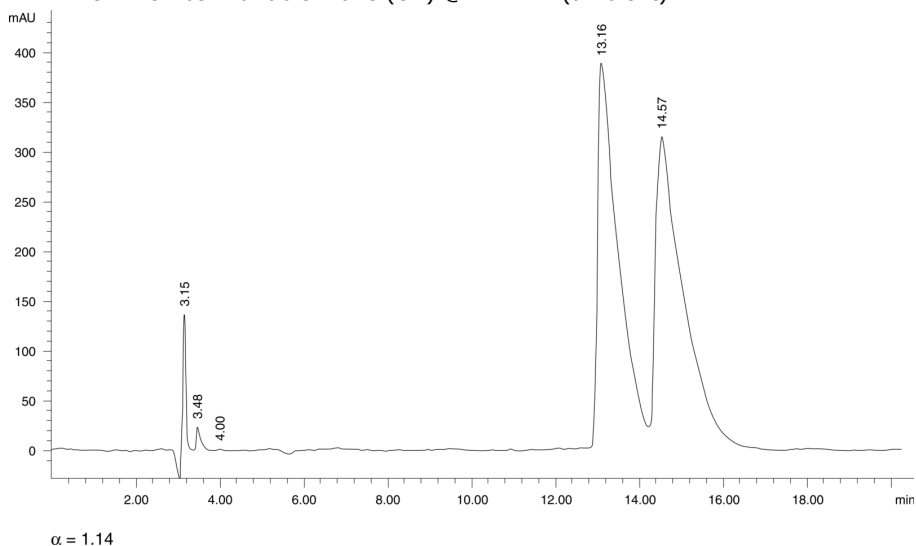


Prilocaine**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (60:35:5) ethanol-TFA 20:1**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 242 nm (ambient)**ANALYTES:****1** Prilocaine

Retention Time: 3.15 min

2 Prilocaine

Retention Time: 13.16 min

3 Prilocaine

Retention Time: 14.57 min

