


Proglumide**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (55:35:10) w/ethanol-TFA (20:1)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 242 nm (ambient) $\alpha = 1.27$ **ANALYTES:****1** Proglumide

Retention Time: 4.64 min

2 Proglumide

Retention Time: 5.16 min

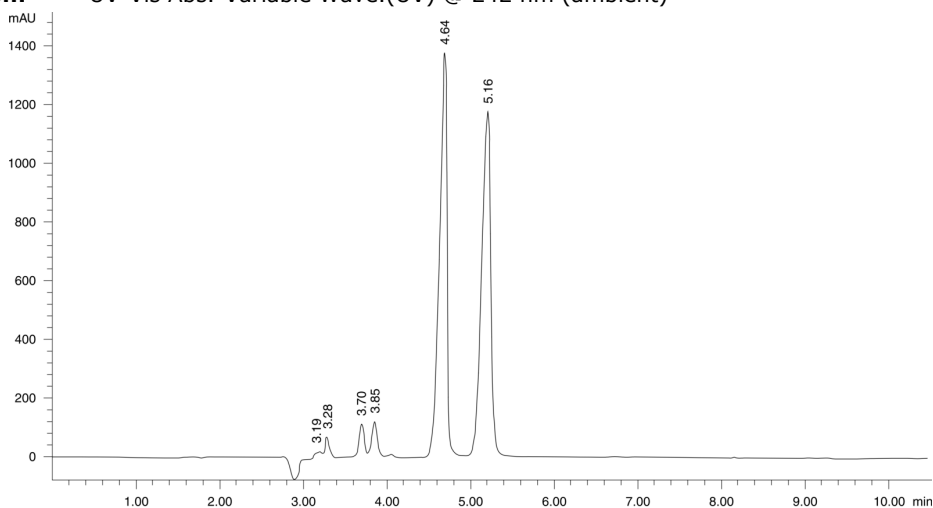