


Promethazine

Column: Chirex® (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3020-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (60:35:5) ethanol-TFA 20:1
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

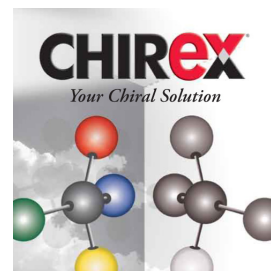

Products used in this application:

CHIREX

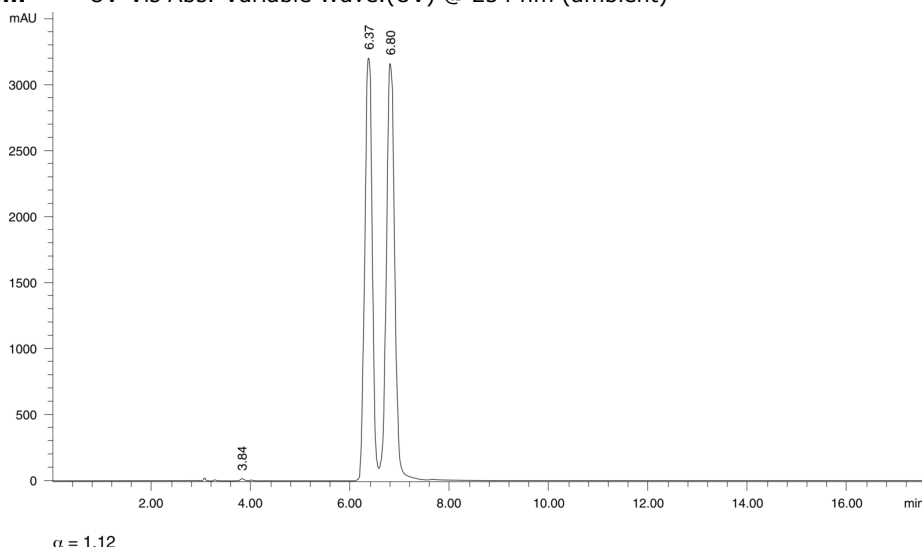

ANALYTES:

- 1 Promethazine
Retention Time: 3.84 min
- 2 Promethazine
Retention Time: 6.37 min
- 3 Promethazine
Retention Time: 6.8 min

