


Promethazine**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (76:20:4) ethanol-TFA 20:1**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**ANALYTES:****1** Promethazine

Retention Time: 32.33 min

2 Promethazine

Retention Time: 34.73 min

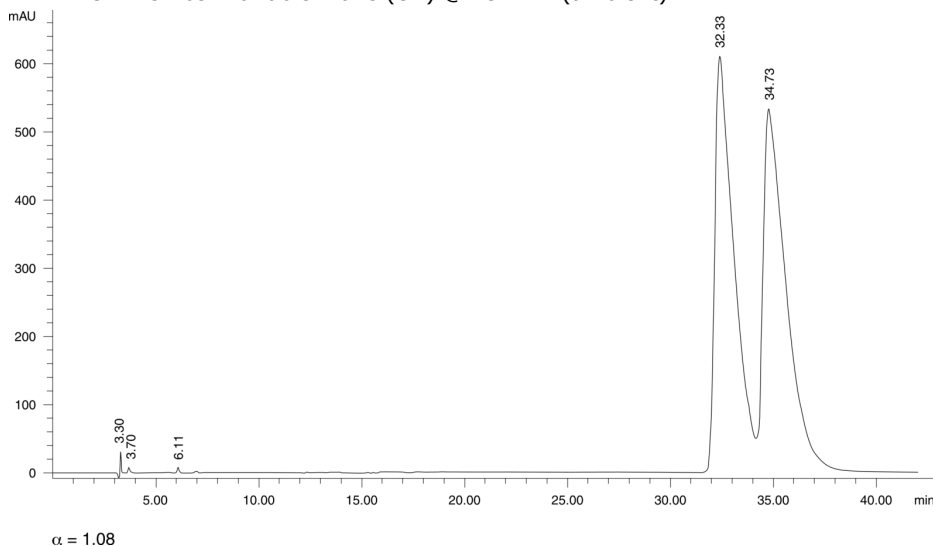