

HPLC Application


ID No.: 13692

Propranolol

Column: Chirex® (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3020-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (76:20:4) ethanol-TFA 20:1
Gradient


Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 290 nm (ambient)

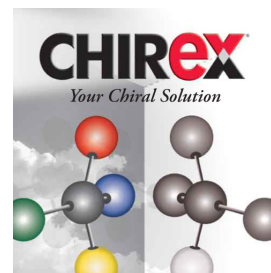

Products used in this application:

CHIREX

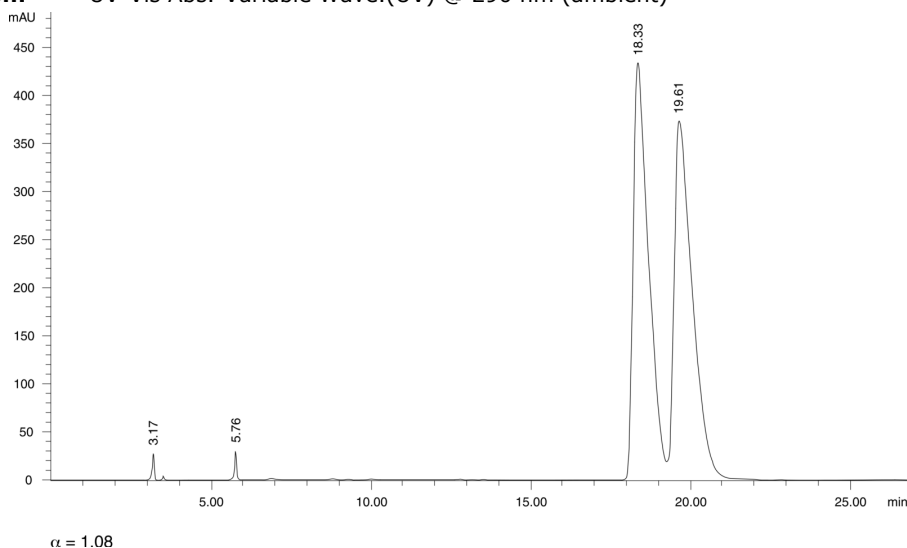

ANALYTES:

- 1 Propranolol
Retention Time: 18.33 min
- 2 Propranolol
Retention Time: 19.61 min