

Terbutaline

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Hexane/Dichloroethane/Ethanol-TFA (55:35:10) ethanol-TFA 20:1


Gradient **Step No.** **Time (min)** **Pct A**

Profile: **1** 0 100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 278 nm (ambient)

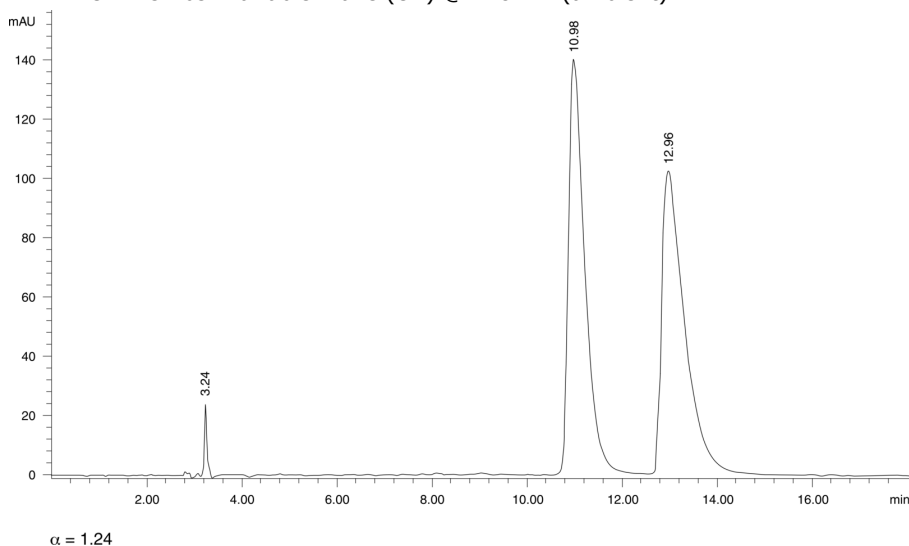

ANALYTES:

1 Terbutaline

Retention Time: 10.98 min

2 Terbutaline

Retention Time: 12.96 min

